Therabody

# SleepMask

User Manual & Unit Warnings

## Table of contents

1.	Product Overview	3
2.	What's in the Box	6
3.	Getting to Know the Device	6
4.	Using the Device	8
5.	How the Device Works	12
6.	Download the Therabody App	12
7.	Care, Cleaning, and Charging	16
8.	Troubleshooting Guide	26

## 1. Product Overview


# SleepMask

Ultra-comfortable 100% blackout eye mask with gentle vibration for better quality sleep.

Better sleep just got easier with Therabody SleepMask device. With an ultra-comfortable 100% blackout design and 3 gentle vibration patterns with adjustable intensity levels, the device helps you fall asleep faster and improve sleep quality, so you wake up feeling refreshed. Designed for all-night comfort in any sleep position, the soft, breathable fabric conforms perfectly to every face shape without any pressure on the eyes, and the adjustable, non-slip strap prevents hair tangling. Plus, with the Therabody app, you can access TheraMind™ sound, an extensive library of music and audio that uses science-backed composition techniques to help you relax and sleep — perfectly paired with the device. Whether you're winding down for the night, traveling, or napping, the device creates the perfect environment for restful sleep anytime, anywhere.

### Intended Use

The Therabody SleepMask device's comfortable, light-blocking design and gentle vibration patterns are intended to help you fall asleep faster, improve quality of sleep, and reduce sleep disturbances so you wake up feeling refreshed and alert.

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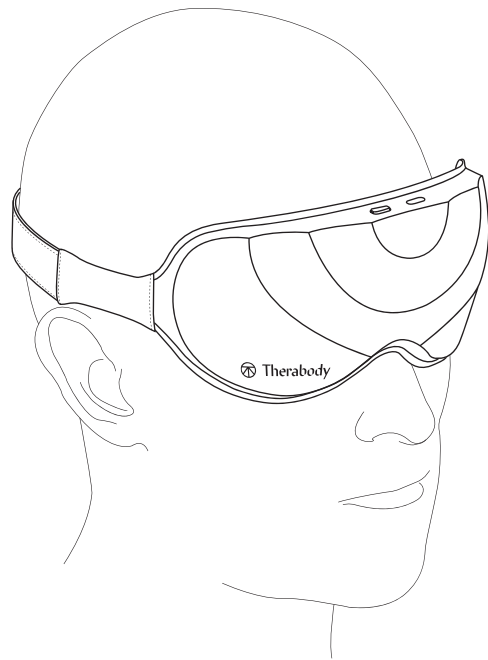
Note: A power adapter is not included with this device. For best performance, use a USB-C power adapter of at least 20W power capacity. Ensure it is from a certified manufacturer and has not suffered any structural damage.
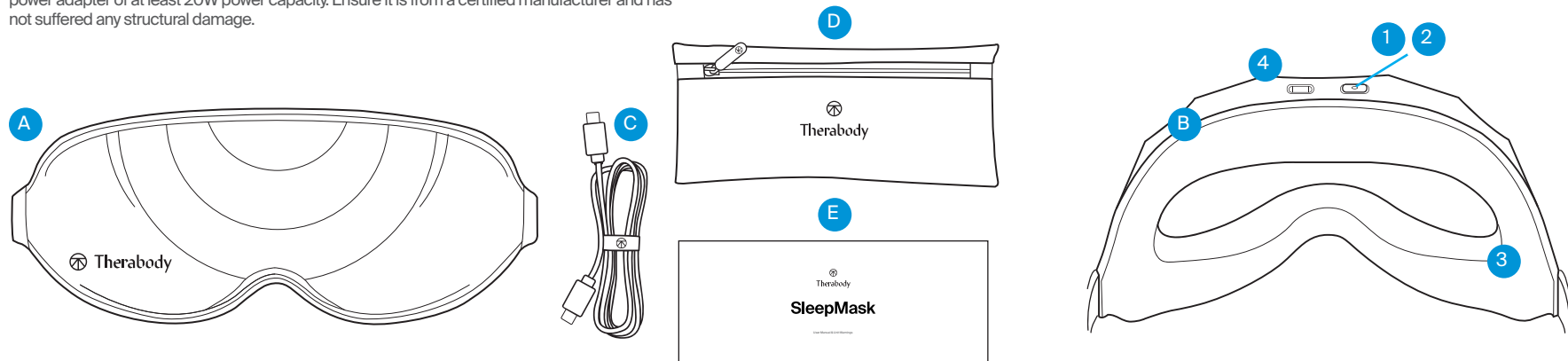
## 3. Getting to Know the Device

1. Multi-functional power/vibration button
2. LED battery status indicator
3. Therabody SleepMask (assembled removable inner and outer control layer)
4. USB-C charging port

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
## 4. Use the Device

4

- 1. Clean the device.** Your device and travel pouch should be cleaned before its first use. Carefully follow the cleaning instructions in the *Cleaning, Care & Charging* section for the inner layer, outer layer and travel pouch to keep your mask in good condition.
- 2. Put on the device.** Once clean, place the device on clean, dry skin comfortably on your face and completely covering your eyes. Remove makeup, eyeglasses, and contact lenses prior to use.
- 3. Adjust the strap.** The strap should sit just above your ears. For a full light block experience, tighten the strap by adjusting the buckle. Moving it to the right makes it tighter and to the left makes it looser. The device should feel snug on your face yet still comfortable.
- 4. Get into your desired sleeping position.** Find a comfortable sleeping position on your back, side or stomach and adjust the device to ensure it is secured.
- 5. Turn on the device.** Firmly press and hold the multi-functional power button located on the top center of the mask. The LED light on the button will turn on when the device is on.
- 6. Select a vibration pattern and intensity.** Press the same button to toggle through the three vibration patterns: Pulse, Wave, and Constant. Each vibration pattern has two intensity levels: low and high, with the first pattern starting at low intensity. Press the button again to switch to high intensity or keep pressing to cycle through each vibration pattern and intensity level. With the stored memory feature, the device will remember the last vibration pattern used so that the next time you turn it on, the same pattern will start.
- 7. Relax and enjoy.** The full session will last for 15 minutes. The LED indicator will turn off 5 seconds after the session starts. The intensity level will gradually decrease as the cycle progresses.
- 8. Turn off the device.** The device will automatically shut off after the session is completed. Press the button at any point during the session to turn off the device. The LED light on the button will turn off when the device is off.
- 9. Lock the device.** When not in use or during transport, lock the device to prevent turning it on accidentally. Quickly press the button three times to initiate the lock function. The device will flash three times when successfully locked. Quickly press three times again to unlock. The device will flash three times when successfully unlocked.

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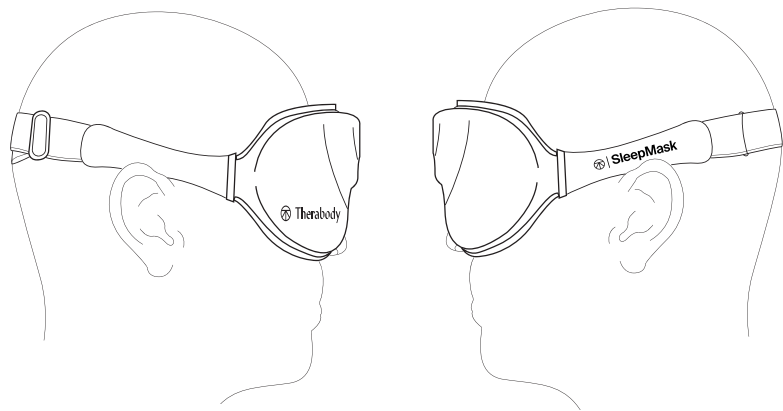## 5. How the Device Works

### 5 6 Modes Overview

The device has three strategically designed vibration patterns, Pulse, Wave and Constant, each with two intensity levels, Low & High (6 total). Vibration can help to promote relaxation and improve sleep quality.

*The Pulse pattern* is a gentle pulsing rhythm like a calm, relaxed heartbeat.

1. Pulse: Low
2. Pulse: High


*The Wave pattern* is a slow rise and fall like a mellow rolling wave, and the inhale and exhale of your breath.

3. Wave: Low
4. Wave: High

*The Constant pattern* is a gentle, steady rhythm — like white noise in physical form.

5. Constant: Low
6. Constant: High

With the stored memory feature, the device will remember the last vibration pattern used so that the next time you turn it on, the same pattern will start.

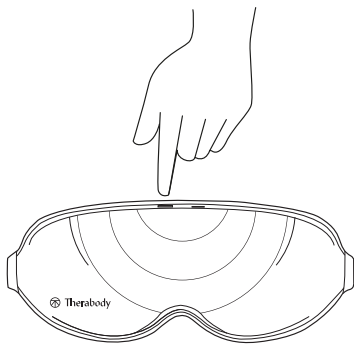*Press and hold the power button to turn on the device. Press the same button to toggle through the vibration patterns and intensity levels.*

## 6. Download the Therabody App

Experience even better sleep with free TheraMind sound therapy content in the Therabody app.

1. Download the Therabody app from the App Store or Play Store.
2. Follow the in-app instructions to create a new app account and verify your email.
3. Once logged in, browse the expansive catalog of TheraMind sound therapy tracks to listen to during your SleepMask treatment. These science-backed audio tracks use special composition techniques to help you relax and sleep.
4. Select a sound that works for you. Whether you prefer the soothing sounds of nature, peaceful piano melodies, or even guided meditations.
5. Save your favorite sounds for easy access.
6. Set a sleep timer so that the sound gently fades as you drift off to sleep.

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


## 7. Care, Cleaning and Charging

### 7 Device Maintenance

The following maintenance instructions are important to ensure that your device continues to work as it was designed. Failure to follow these instructions may cause your device to stop working.

#### Assemble and Disassemble the Device

1. Lay the device flat, with the inside facing up towards you.
2. Move the elastic strap to the front of the device, so it is out of the way.
3. Hold one end of the eye mask's removable inner layer and gently peel it off to detach it from the outer control layer.
4. To reattach, hold both ends of the removable eye mask to reattach it to the front layer. There is a magnet in the center to ensure the layers are correctly attached.

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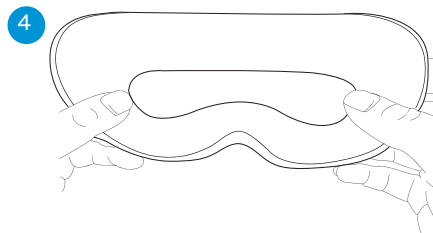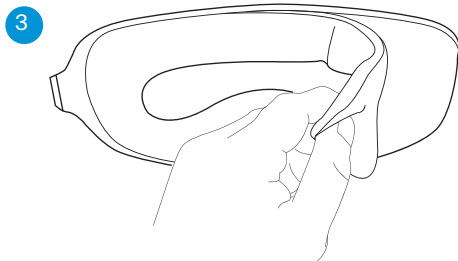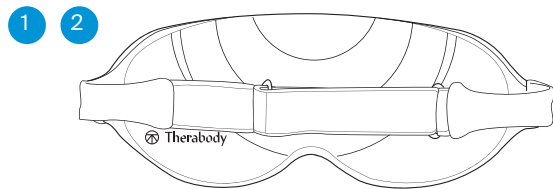## Cleaning the Device

Your device and travel pouch should be cleaned before its first use. For ongoing care, clean at least once a week, similar to your bedding. Carefully follow the cleaning instructions below to ensure longevity of the device and its components.

### 7 Removable Inner layer + Travel pouch

1. Carefully remove the inner layer from the device following the instructions above. Place **ONLY** the removable inner layer into the washing machine.
2. Machine wash the removable inner layer on cold with mild detergent using the delicate wash cycle.
3. Once the wash cycle is complete, remove the inner layer from the washing machine and lay it flat to air dry completely in a cool and dry area (this can take up to 8 hours or more, depending on conditions). **Wait until completely dry before reassembling and using the device.** DO NOT place the inner layer in the dryer, as the heat will damage the fabric. DO NOT wring, twist, or squeeze the inner layer to dry.
4. Follow the same cleaning instructions above to wash the travel pouch. Hand-wash after 75 cycles.

### Outer Control Layer

1. The outer control layer of the device is **NOT** waterproof, hand-washable, or machine-washable. Do not allow excess liquids or water to enter the device.
2. To clean the outer control layer, wipe the fabric with a non-abrasive, lint-free cloth.
3. Lightly dampen a non-abrasive, lint-free cloth with water or an alcohol-based solution to clean if necessary. Do not allow water to touch or enter the power button or USB-C port.
4. Dry with a non-abrasive, lint-free cloth.

\*Note: For the best experience, we recommend replacing the removable inner layer annually, available at [www.therabody.com](http://www.therabody.com).

## Charging the Device

1. The device is USB-C enabled and comes with a charging cable. Before charging the device, ensure that the device is powered off.
2. Connect one end of the charging cable to the USB-C charging port on the device. It is located in the top center next to the power button. Connect the other end of the charging cable to a power adapter.
3. The LED light will indicate the battery status when the device is on and charging. Press the button to check the battery status while the device is off.
  - Flashing orange: Low battery
  - Flashing blue: Half battery
  - Flashing green: High battery
  - Solid green: Full battery

\*\*Note: If using a charger from a third party, ensure that it is from a certified manufacturer and has not suffered any structural damage. Do not try to use the device while charging.

Label	Description	Location of Label
	Cold wash this layer, air dry only	Removable Inner layer + Travel pouch
	Hand Wash	Removable Inner layer + Travel pouch
	Do not dry-clean	Removable Inner layer + Travel pouch
	Do not tumble dry	Removable Inner layer + Travel pouch
	Do not bleach	Removable Inner layer + Travel pouch
	Do not iron	Removable Inner layer + Travel pouch
	Max washing temperature 20°C normal process	Removable Inner layer + Travel pouch
	Do not wash this layer	Outer Control Layer

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## 8. Troubleshooting

SCENARIO	SYMPTOM	POTENTIAL CAUSES	POTENTIAL SOLUTIONS
<p><b>Device does not turn on.</b></p>	<p>Battery indicator LED doesn't turn on</p>	<p>Battery ran out or No battery power in the device</p>	<p>Connect the USB-C cable (included in box) to the device and plug in the other end to a power outlet. Verify that the device is charging by checking proper connection between the power adapter, the USB-C cable, and the USB-C port on the device.</p>
		<p>USB-C cable is damaged</p>	<p>Ensure you are using the USB-C to USB-C cable provided, or one from a certified manufacturer.</p>
		<p>Internal malfunction</p>	<p>If the device is not charging after verifying all the above, Contact Therabody for further assistance. Visit <a href="https://www.therabody.com/us/en-us/support/support.html">https://www.therabody.com/us/en-us/support/support.html</a> for the available contact methods.</p>
	<p>Power LED can light on but the device not work</p>	<p>Battery ran out or No battery power in the device</p>	<p>Charge the device and try again</p>

SCENARIO	SYMPTOM	POTENTIAL CAUSES	POTENTIAL SOLUTIONS
Slow charging	Slow charging or incompletely charging	Adapter issue or Adapter power too low	Ensure you are using the USB-C to USB-C cable provided, or one from a certified manufacturer.
			If the issue still persists, Contact Therabody for further assistance. Visit <a href="https://www.therabody.com/us/en-us/support/support.html">https://www.therabody.com/us/en-us/support/support.html</a> for the available contact methods.
Device does not turn on	Device does not turn on and the LED flashes white three times	Device is locked	Press the button three times. The device will vibrate twice and the LED will flash green twice. Then, the device will unlock.
			After verifying the above, if the device does not turn on, the device may be damaged. Contact Therabody for further assistance. Visit <a href="https://www.therabody.com/us/en-us/support/support.html">https://www.therabody.com/us/en-us/support/support.html</a> for the available contact methods.
Vibration function is not working	Vibration does not turn on even after vibration button is pressed	Internal Malfunction	If the vibration does not turn on after verifying the above, Contact Therabody for further assistance. Visit <a href="https://www.therabody.com/us/en-us/support/support.html">https://www.therabody.com/us/en-us/support/support.html</a> for the available contact methods.
Mechanical Noise	Rattling noise when powered ON	Internal Malfunction	Contact Therabody for further assistance. Visit <a href="https://www.therabody.com/us/en-us/support/support.html">https://www.therabody.com/us/en-us/support/support.html</a> for the available contact methods.

DO NOT attempt to troubleshoot by disassembling the device as this will void the warranty and may cause injury.

**Store and Transport the Device**

To store and transport the device, lay it flat in its travel pouch to not damage the internal components.

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# Warnings and Guidance (Precautions and Contraindications)

## Background

Therabody products are designed to unlock the body's natural ability to achieve health and well-being. Through science and technology, the Therabody portfolio of products allows people to access the benefits of different natural phenomena to meet both their needs and preferences. There will be times when it is advisable to modify how devices are used (precautions) or times when it is not appropriate to use certain devices (contraindications). Read the following safety information for the device in its entirety prior to use.

## Important Safety Information

### General device Use

This device is intended for use by people in good health. This device is contraindicated against and should not be used by or on anyone with a history of epilepsy, seizures or cardiopathy.

The the device is not recommended for anyone with an electronic implanted device (such as a pacemaker), cardiac arrhythmia, tumors, or acute episodes of inflammatory diseases. The device is not recommended for those who have arteriosclerosis, thromboses, or implants in the body region being treated.

Please consult your physician prior to using the device if you are pregnant and/or nursing. **Immediately stop using the device at the first sign of discomfort.** The device IS NOT a toy. If you have any medical considerations please consult your doctor before using the device.

## Safety, Precautions, and Contraindications

### Specific Therapies

These recommendations are derived from consultation with medical experts and published research regarding precautions and contraindications as of the printing date. For up-to-date information, please visit us online at <https://www.therabody.com/us/en-us/precautions-and-contraindications.html>.

## Precautions:

Due care is required in these circumstances and device use may need to be modified. Consult with a medical professional if you currently have or suspect you may have any of the following conditions or if you have any questions.

- Recent injury or surgery
- Hypertension (controlled)
- Abnormal sensations (e.g., numbness)
- Sensitivity to pressure
- Bony prominences or regions
- Mild peripheral neuropathy
- It's recommended that you remove contact lenses before using the device
- This device should fit tightly but should NOT restrict circulation in any way

## Contraindications:

The following are circumstances where the potential risks may outweigh the benefits. Consult a medical professional before use.

- Skin rash, open wounds, blisters, local tissue inflammation, infections, bruises, or tumors
- Hypertension (uncontrolled)
- Vascular/arterial and circulation diseases, including blood clots, phlebitis, thromboses, and varicose veins
- Acute or severe cardiac, liver, or kidney disease
- Severe peripheral neuropathy or other causes of sensory impairment
- Bleeding disorders
- Connective tissue disorders
- Medications that may alter or thin the blood or alter sensations
- Direct pressure over surgical site or hardware
- Pacemaker or ICD
- History of embolism

- Allergy to device material (polyester, spandex, PU foam)
- Do not use in the presence of unexplained pain

#### Limited Warranty and Customer Service

For full warranty information, please visit [www.therabody.com/warranty](http://www.therabody.com/warranty).

Customers who are in need of product support should visit

<https://www.therabody.com/us/en-us/support/support.html> for the available contact methods.

#### FDA-Registered

Limited Warranty Only With Authorized Retailer Purchase

© 2024 Therabody, Inc. All Rights Reserved.

Patents at [www.therabody.com/patents](http://www.therabody.com/patents)

#### Manufactured and Distributed by:

Therabody, Inc.

1640 S Sepulveda, Suite 300

Los Angeles, CA 90025

In the Australian market, the Therabody SleepMask warranty is 1 year.

For Australian Market Only: The purchase of this device is accompanied by a limited warranty. For complete details about the rights afforded to you under this warranty, visit <https://therabody.com/anz/en-au/warranty-support/>.

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## UNIT WARNINGS

READ ALL INSTRUCTIONS AND CAUTIONARY MARKINGS IN THIS MANUAL, ON THE CHARGER, AND ON THE DEVICE BEFORE USING OR CHARGING THE DEVICE AND ITS ACCESSORIES.

When using the device, the following precautions should always be adhered to:

1. **USE ONLY AS INSTRUCTED.** Use only Therabody recommended accessories and replacement parts. This device is not intended to diagnose, mitigate, or treat diseases. The device is an advanced mechanical tool with electric components. If the device and its accessories are not used or maintained properly, there is a risk of fire, electric shock, or injury. Failure to use and maintain the device and corresponding attachments in accordance with the instructions outlined in this manual will void the warranty and may result in product damage or physical injury. The device is intended for commercial and home use.
2. **NOT FOR CHILDREN.** The device is not intended for use by young children or persons with reduced physical, sensory, or reasoning capabilities or lack of experience and knowledge, unless given supervision or instruction concerning use of the device by a person responsible for their safety. The device is not a toy. Do not play with, bend, or pull the electrical components. Advise children not to play with the device or accessories. Adult supervision should be provided for those under 18 using the device.
3. **CHARGING.** If your device doesn't turn on or the battery indicator displays a low battery level, please charge before first use. Charge the battery using the USB-C cable (included). When you use the charging cable (included), make sure its USB-C connector is fully inserted into a compatible power adapter before you plug the adapter into a power outlet. You can also charge the battery with third party cables that are compliant with USB-C and with applicable country regulations and international and regional safety standards. Other adapters may not meet applicable safety standards, and charging with such adapters could pose a risk of death or injury. Do not force the cable into place. Fully charge the device at least once every six months to prevent device and battery damage. Unplug the power adapter and charging cable when not in use.
4. **DO NOT OVERCHARGE.** Do not leave the device connected to the power adapter for more than one hour after the battery has been fully charged. The battery includes a protection system to avoid the risk of overcharging. However, overcharging may reduce its life over time.
5. **CHARGING LOCATIONS.** Ensure that the battery status LED is OFF before charging. Insert the USB-C cable into the port (behind the power button) to charge the device. Do not use the device while charging. Avoid skin contact when the USB-C is plugged into a power source as it may cause discomfort or injury. Do not sleep or sit on the USB-C connector. The device should be charged indoors in a well-ventilated, dry location. Do not charge the device outdoors, in a bathroom, or within 10 feet (3.1 meters) of a bathtub or pool. Do not use the device or charger on wet surfaces, and do not expose the charger to moisture, rain, or snow. Do not use the device or its compatible charger in the presence of explosive atmospheres (gaseous fumes, dust, grain, metal powders, or flammable materials). Sparks may be generated, possibly causing a fire.
6. **UNPLUGGING THE CHARGER.** Pull the plug, not the cord, to reduce the risk of damage to the electrical plug and cord. Never carry the charger by its cord. Keep the cord away from heat, oil, and sharp edges. Do not stretch the charger cord or place the cord under strain. Do not handle the charger, terminals, or the device with wet hands. Always unplug this appliance from the electrical outlet immediately after using and before cleaning. Store cable to ensure it is not stepped on, tripped over, or otherwise subjected to damage or stress. Do not use a power adapter or charging cable that has received a sharp blow, been dropped, run over, or damaged

in any way. For long-term storage, store with a fully charged battery.

7. **USAGE AND SKIN SENSITIVITIES.** Adjust the strap to fit comfortably on your head. The device should perfectly align with your temples and sit comfortably on your nose bridge. Always adjust the device while the device is OFF. Long press the Power Button to turn off the device. Make sure the Battery status LED is OFF. Discontinue the use of the device if you feel any discomfort. Use the device on clean, dry skin. Remove makeup before device use. Some people may experience reactions to certain materials used in wearable items that are in prolonged contact with their skin. This can be due to allergies, environmental factors, extended exposure to irritants like sweat or other causes. You may be more likely to experience irritation from any wearable device if you have allergies or other sensitivities. If you have known sensitivities, take special care when using the device and accessories. If you experience redness, swelling, itchiness, or any other irritation or discomfort on your skin around, or beneath, your device, remove the device and consult your physician before resuming wear. Continued use, even after symptoms subside, may result in renewed or increased irritation.

8. **DO NOT OPERATE UNDER BLANKET AND PILLOW.** Excessive heating can occur and cause fire, electric shock, or injury.

9. **DEVICE CARE AND SERVICE.** Remove eyeglasses before using the device. It's recommended that you remove contact lenses before using the device. Handle the device with care. The device contains electronic components that can be damaged if dropped, burned, punctured, or crushed. Do not disassemble the device or attempt to repair the device yourself. Disassembling the device may damage it, result in loss of water resistance, and may cause injury. If the device is damaged or if malfunctions occur, contact Therabody Customer service. The warranty will be void if the device, batteries, or charger are disassembled or if any parts have been removed. Do not use if damaged as these may cause injury. Avoid heavy exposure to dust or sand. Do not puncture or damage the device. Puncturing the device can lead to incorrect operation.

10. **CLEANING INSTRUCTIONS.** The outer electronic layer is NOT waterproof or machine washable. Clean the outer layer, elastic strap, and side fabric with a lint-free damp cloth to remove stains, dirt, makeup, ink, acids, or acidic food. Do not use a cloth soaked in water to clean the device. Only the inner layer is machine washable. Carefully remove the inner layer from the device and place ONLY the inner layer into the washing machine. Machine wash the inner layer on cold with mild detergent using the delicate wash cycle. Once the wash cycle is complete, remove the inner layer from the washing machine and lay it flat to air dry completely (this can take up to 8 hours or more depending on conditions). DO NOT place the inner layer in the dryer, as the heat may damage the fabric. Wait until completely dry before reassembling the device. DO NOT wring, twist, or squeeze the inner layer to dry. A properly cleaned device should have no visible signs of debris or moisture.

11. **STORING THE DEVICE AND ACCESSORIES.** Store in a cool, dry place. Only charge the device when the ambient temperature is between 40°C/104°F and 10°C/50°F. Do not store the device or charging cable where temperatures may exceed 40°C/104°F, such as in direct sunlight, in a vehicle, or in a metal building during the summer.

12. **DEVICE DISPOSAL.** This device contains a lithium-ion battery, and care must be taken upon disposal of the device. Before disposal of this device, please review your local laws and requirements surrounding Lithium-Ion Battery disposal. The preferred method of disposal is recycling the whole device.

### Reporting adverse events to FDA

MedWatch is the Food and Drug Administration's (FDA) program for reporting serious reactions, product quality problems, therapeutic inequivalence/failure, and product use errors with human medical products, including drugs, biologic products, medical devices,

dietary supplements, infant formula, and cosmetics.

If you think you or someone in your family has experienced a serious reaction to a medical product, you are encouraged to take the reporting form to your doctor. Your health care provider can provide clinical information based on your medical record that can help the FDA evaluate your report. However, we understand that for a variety of reasons, you may not wish to have the form filled out by your healthcare provider, or your health care provider may choose not to complete the form. Your health care provider is not required to report to the FDA. In these situations, you may complete the Online Reporting Form yourself.

You will receive an acknowledgement from the FDA when your report is received. Reports are reviewed by FDA staff. You will be personally contacted only if we need additional information.

### Submitting adverse event reports to FDA

Use one of the methods below to submit voluntary adverse event reports to the FDA:

1. Report Online at: [www.accessdata.fda.gov/scripts/medwatch/index.cfm?action=reporting.home](http://www.accessdata.fda.gov/scripts/medwatch/index.cfm?action=reporting.home)
2. Consumer Reporting Form FDA 3500B. Follow the instructions on the form to either fax or mail it in for submission. For help filling out the form, see MedWatchLearn. The form is available at: [www.fda.gov/downloads/aboutFDA/reportsmanualsforms/forms/ucm349464.pdf](http://www.fda.gov/downloads/aboutFDA/reportsmanualsforms/forms/ucm349464.pdf)
3. Call FDA at 1-800-FDA-1088 to report by telephone.
4. Reporting Form FDA 3500 commonly used by health professionals. The form is available at: [www.fda.gov/downloads/aboutFDA/reportmanualsforms/forms/ucm163919.pdf](http://www.fda.gov/downloads/aboutFDA/reportmanualsforms/forms/ucm163919.pdf)

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment

Label	Description	Location of Label
	Do not dry-clean	Removable Inner layer + Travel pouch
	Do not tumble dry	Removable Inner layer + Travel pouch
	Do not bleach	Removable Inner layer + Travel pouch
	Do not iron	Removable Inner layer + Travel pouch
IP 22	Degree of protection against ingress of water	On rating label
	Read instructions before use	On rating label
	Level of protection type BF applied part	On rating label
	Therabody, Inc. 1640 S. Sepulveda Blvd. Suite 300 Los Angeles, CA 90025	On rating label

Label	Description	Location of Label
	Unique Device Identification (UDI)	Packaging
	Separate collection for waste electrical and electronic equipment. Note: For more information about disposal of equipment, its parts and accessories, please contact your local distributor.	On rating label
	UKCA Mark for UK	On garment label
	IC (Canada)	On garment label
	RCM	On garment label

This device is type BF applied equipment . It complies with Medical Electrical Safety Standards (IE60601-1). This device also complies with Medical EMC Standard (IEC 60601-1-2).The has been tested and found to comply with the electromagnetic compatibility (EMC) limits for medical devices to IEC 60601-1-2: 2014. These limits are designed to provide reasonable protection against harmful interference in atypical medical installation.

**WARNING:**

1\* **WARNING:** Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

2\* **WARNING:** Use of accessories, transducers and cables other than those specified or provided by the manufacturer this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in proper operation.

3\* **WARNING:** Portable RF communications equipment (including peripherals such as antenna cables) should be used no closer than 30 cm (12 inches) to any part of the ME EQUIPMENT, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result

Guidance and manufacturer's declaration – electromagnetic emission

The Sleep Mask is intended for use in the electromagnetic environment specified below. The customer or the user of Sleep Mask should assure that it is used in such an environment.

Emissions test	Compliance	Electromagnetic environment - guidance
RF emissions CISPR 11	Group 1	The Sleep Mask uses RF energy only for its internal function. There for, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The Sleep Mask suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations flicker emissions IEC 61000-3-3	Not applicable	

## Guidance and manufacturer's declaration – electromagnetic immunity – for all EQUIPMENT and SYSTEMS

The Sleep Mask is intended for use in the electromagnetic environment specified below. The customer or the user of the Sleep Mask should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment -guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrostatic transient / burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	N/A	Mains power quality should be that of a typical commercial or home environment.
Surge IEC 61000-4-5	± 1 kV differential mode ± 2 kV common mode	N/A	Mains power quality should be that of a typical commercial or home environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT; 0,5 cycle g) At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0 % UT; 1 cycle and 70 % UT; 25/30 cycles Single phase: at 0° 0 % UT; 250/300 cycle	N/A	Mains power quality should be that of a typical commercial or home environment. If the user of the Sleep Mask requires continued operation during power mains interruptions, it is recommended that the Sleep Mask be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or home environment.

NOTE      UT is the a. c. mains voltage prior to application of the test level.

## Guidance and manufacturer's declaration – electromagnetic immunity - for EQUIPMENT and SYSTEM

The Sleep Mask is intended for use in the electromagnetic environment specified below. The customer or the user of the Sleep Mask should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment –guidance
<p>Conducted RF IEC 61000-4-6</p>	<p>3 Vrms 150 kHz to 80 MHz</p> <p>6 V in ISM and amateur radio bands between 0,15 MHz and 80 MHz</p>	N/A	<p>Portable and mobile RF communications equipment should be used no closer to any part of the Sleep Mask, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended separation distance</p> $d = \left[ \frac{3.5}{V_1} \right] \sqrt{P}$ $d = \left[ \frac{12}{V_2} \right] \sqrt{P}$
<p>Radiated RF IEC 61000-4-3</p>	<p>10 V/m 80 MHz to 2.7 GHz</p> <p>385MHz-5785MHz Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communication equipment (Refer to table 9 of IEC 60601-1-2:2014)</p>	<p>10 V/m 80 MHz to 2.7 GHz</p> <p>385MHz-5785MHz Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communication equipment (Refer to table 9 of IEC 60601-1-2:2014)</p>	$d = \left[ \frac{3.5}{E_1} \right] \sqrt{P} \quad 80 \text{ MHz to } 800 \text{ MHz}$ $d = \left[ \frac{7}{E_1} \right] \sqrt{P} \quad 800 \text{ MHz to } 2.7 \text{ GHz}$ <p>where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m).b</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,a should be less than the compliance level in each frequency range.b</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol: </p>

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic is affected by absorption and reflection from structures, objects and people.

a The ISM (industrial, scientific and medical) bands between 150 kHz and 80 MHz are 6,765 MHz to 6,795 MHz; 13,553 MHz to 13,567 MHz; 26,957 MHz to 27,283 MHz; and 40,66 MHz to 40,70 MHz. The amateur radio bands between 0,15 MHz and 80 MHz are 1,8 MHz to 2,0 MHz, 3,5 MHz to 4,0 MHz, 5,3 MHz to 5,4 MHz, 7 MHz to 7,3 MHz, 10,1 MHz to 10,15 MHz, 14 MHz to 14,2 MHz, 18,07 MHz to 18,17 MHz, 21,0 MHz to 21,4 MHz, 24,89 MHz to 24,99 MHz, 28,0 MHz to 29,7 MHz and 50,0 MHz to 54,0 MHz.

b Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Sleep Mask is used exceeds the applicable RF compliance level above, the Sleep Mask should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the Sleep Mask.

c Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3V/m.

## Recommended separation distances between portable and mobile RF communications equipment and the EQUIPMENT or SYSTEM - for EQUIPMENT and SYSTEMS

The Sleep Mask is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Sleep Mask can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Sleep Mask as recommended below, according to the maximum output power of the communications equipment

### Separation distance according to frequency of transmitter - m

Rated maximum output of transmitter W	150 kHz to 80 MHz outside ISM and amateur radio bands	150 kHz to 80 MHz in ISM and amateur radio bands	80 MHz to 800 MHz	800 MHz to 2.7 GHz
	$d = \left[ \frac{3.5}{V_1} \right] \sqrt{P}$	$d = \left[ \frac{12}{V_2} \right] \sqrt{P}$	$d = \left[ \frac{3.5}{E_1} \right] \sqrt{P}$	$d = \left[ \frac{7}{E_1} \right] \sqrt{P}$
0.01	0.12	0.20	0.035	0.07
0.1	0.38	0.63	0.11	0.22
1	1.2	2.00	0.35	0.70
10	3.8	6.32	1.10	2.21
100	12	20.00	35	70

For transmitters rated at a maximum output power not listed above the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

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#### ISED compliance statement Wireless requirement

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

(1) This device may not cause interference.

(2) This device must accept any interference, including interference that may cause undesired operation of the device.

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